

Shijiazhuang Hongray Group Co., Ltd.

South Tongda Rd., East Dist., Jinzhou City, Hebei, 052260, China

The following merchandise was submitted and identified by the applicant as:

Product Description: Disposable Nitrile Glove (blue beaded ambidextrous)

Style/Item No.: NPF2005-XL

We have tested the submitted sample(s) as requested and the following results were obtained:

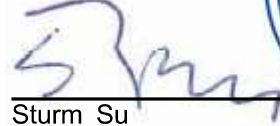
Test Requested: Viral Penetration Test (Refer to ASTM F 1671-2013)

Test Method & Result: ---See following sheet(s)---

Date of Receipt : MAY 06, 2019

Testing Period : MAY 06, 2019 ~ MAY 24, 2019

Signed for and on behalf of
SGS Taiwan Ltd.


Sturm Su
Asst. Manager



This document is issued by the Company subject to its General Conditions of Service printed overleaf, available on request or accessible at <http://www.sgs.com/en/Terms-and-Conditions.aspx> and, for electronic format documents, subject to Terms and Conditions for Electronic Documents at <http://www.sgs.com/en/Terms-and-Conditions/Terms-e-Documents.aspx>. Attention is drawn to the limitation of liability, indemnification and jurisdiction issues defined therein. Any holder of this document is advised that information contained hereon reflects the Company's findings at the time of its intervention only and within the limits of client's instruction, if any. The Company's sole responsibility is to its Client and this document does not exonerate parties to a transaction from exercising all their rights and obligations under the transaction documents. This document cannot be reproduced, except in full, without prior written approval of the Company. Any unauthorized alteration, forgery or falsification of the content or appearance of this document is unlawful and offenders may be prosecuted to the fullest extent of the law. Unless otherwise stated the results shown in this test report refer only to the sample(s) tested

Test Method & Result:

Test Items		Results				Test Method
		Step			Assay	
		1	2	3	(PFU/mL)	
Penetration of Phi-X174 Bacteriophage	Control (+)	None Seen	Seen	--	--	ASTM F1671-2013
	Control (-)	None Seen	None Seen	None Seen	<1	
	Specimen 1	None Seen	None Seen	None Seen	<1	
	Specimen 2	None Seen	None Seen	None Seen	<1	
	Specimen 3	None Seen	None Seen	None Seen	<1	

Note:

- Testing procedure of Penetration of Phi-X174 Bacteriophage:
 - Step 1: Observe for 5 min at 0 psi.
 - Step 2: Slowly increase the pressure to 2.0 psi at the rate of no more than 0.5 psi/s, keep the pressure at 2.0 psi, observe for 1 min.
 - Step 3: Slowly decrease the pressure to 0 psi at the rate of no more than 0.5 psi/s, observe for 54 min.
- 1 psi = 6.895 kPa
- PFU: Plaques Forming Unit
- 2.0 E+2 is equal to 200, 1.3 E+4 is equal to 13000, etc.
- A retaining screen is used.
- Sample description is given by the client: Disposable Nitrile Glove.

Remark:

- The test report is the test results issued by the testing institution as requested by the consignor, it shall not determine the legitimacy of the product.
- This testing was performed by Taiwan Textile Research Institute.

TEST REPORT

– Picture(s) –



Photo A: Appearance of the sample(s)

---End of Report---